



EMERALD
S C I E N T I F I C

Foxx
Life Sciences
Our Innovation. Your Opportunity.

PUREGRIP® Glass Bottles



- Available in clear and amber options
- Sizes range from 100 mL - 20 L
- Highest Polypropylene cap temperature rating in the industry, 142°C
- Industry standard GL45 thread, compatible with Autofil™ Bottle-Top Filters and EZBio® Single-Use Cap Assemblies
- Large, visible writing surface on bottle for additional labeling; white graduations for identification of approximate volume; retrace code for identification
- BOROSIL® 3.3 Borosilicate Glass is mechanically strong, chemically resistant, and withstands thermal shock
- Bottle, VersaCap®, and Pouring Ring are autoclavable and complies with USP-VI


PUREGRIP® Round Glass Media Storage Bottles combine the leading quality of BOROSIL® 3.3 Borosilicate Glass with the patented advantages of the new Foxx GL45 VersaCap® to produce a true leap forward in laboratory glassware. Ideal for mixing, storage, and transportation of media, buffers, chemicals, and solvents, **PUREGRIP® Bottles** offer versatility across a wide variety of laboratory needs. And with Foxx's revolutionary GL45-threaded polypropylene VersaCap®, **PUREGRIP® Bottles** are substantially easier to hold securely, improving lab safety and user experience.

PUREGRIP® Bottles are autoclavable, reusable, and graduated in mL. A large writing space allows for visible labeling, and an extra frosted surface on the VersaCap® offers a convenient second location for writing. Manufactured from top-quality glass, **PUREGRIP® Bottles** are highly resistant to thermal shock. And with an industry-leading maximum temperature rating of 142°C, the VersaCap® stands up to even the toughest conditions.

- 1 Frosted surface on top of cap for convenient labeling
- 2 Patented SECUREGRASP technology built into cap for ultimate ease of handling, especially with gloved hand
- 3 USP Class VI, Animal-Free Polypropylene VersaCap® and Pouring Ring for leakproof, drip-free experience
- 4 Sloped shoulders for smooth, natural pouring motion
- 5 Glass bead indicates full capacity level
- 6 Large, visible writing surface on bottle for additional labeling; white graduations for identification of approximate volume; retrace code for identification
- 7 Easy to read graduation marks



PUREGRIP® is a registered trademark of Foxx Life Sciences.

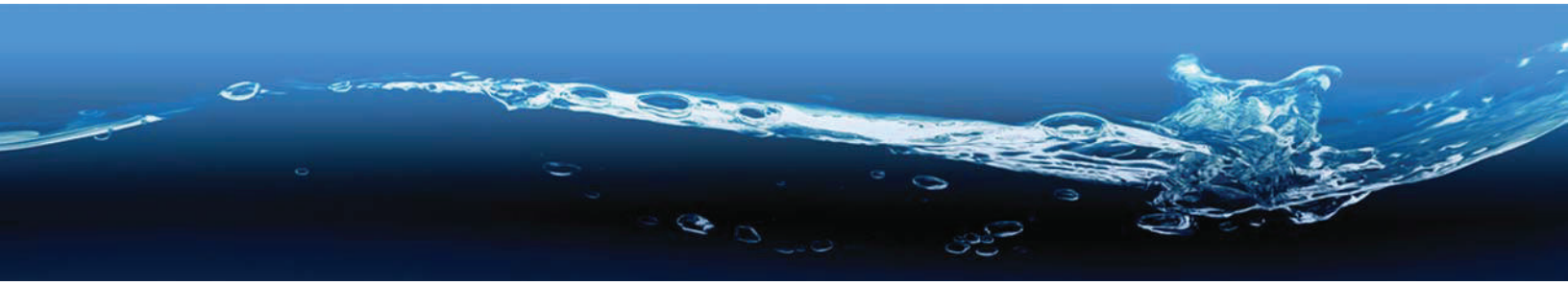

| | Product Number | Description | Total Capacity (mL) | Approx. Dimensions (Dia. x H) | Thread Type | Quantity Per Case |
|-------|----------------|---|---------------------|-------------------------------|-------------|-------------------|
| Clear | 1501016-FLS | PUREGRIP® Reusable Glass Media Storage Bottles with VersaCap® | 100 | 56 x 103 | GL45 | 10 |
| | 1501021-FLS | PUREGRIP® Reusable Glass Media Storage Bottles with VersaCap® | 250 | 70 x 141 | GL45 | 10 |
| | 1501024-FLS | PUREGRIP® Reusable Glass Media Storage Bottles with VersaCap® | 500 | 86 x 179 | GL45 | 10 |
| | 1501029-FLS | PUREGRIP® Reusable Glass Media Storage Bottles with VersaCap® | 1,000 | 101 x 229 | GL45 | 10 |
| | 1501030-FLS | PUREGRIP® Reusable Glass Media Storage Bottles with VersaCap® | 2,000 | 136 x 265 | GL45 | 10 |
| | 1501031-FLS | PUREGRIP® Reusable Glass Media Storage Bottles with VersaCap® | 3,000 | 158 x 298 | GL45 | 2 |
| | 1501033-FLS | PUREGRIP® Reusable Glass Media Storage Bottles with VersaCap® | 5,000 | 183 x 326 | GL45 | 1 |
| | 1501038-FLS | PUREGRIP® Reusable Glass Media Storage Bottles with VersaCap® | 10,000 | 235 x 408 | GL45 | 1 |
| | 1501040-FLS | PUREGRIP® Reusable Glass Media Storage Bottles with VersaCap® | 20,000 | 296 x 503 | GL45 | 1 |
| Amber | 1519016-FLS | PUREGRIP® Reusable Amber Glass Media Storage Bottles with VersaCap® | 100 | 56 x 103 | GL45 | 10 |
| | 1519021-FLS | PUREGRIP® Reusable Amber Glass Media Storage Bottles with VersaCap® | 250 | 70 x 141 | GL45 | 10 |
| | 1519024-FLS | PUREGRIP® Reusable Amber Glass Media Storage Bottles with VersaCap® | 500 | 86 x 179 | GL45 | 10 |
| | 1519029-FLS | PUREGRIP® Reusable Amber Glass Media Storage Bottles with VersaCap® | 1,000 | 101 x 229 | GL45 | 10 |
| | 1519030-FLS | PUREGRIP® Reusable Amber Glass Media Storage Bottles with VersaCap® | 2,000 | 136 x 265 | GL45 | 10 |
| | 1519031-FLS | PUREGRIP® Reusable Amber Glass Media Storage Bottles with VersaCap® | 3,000 | 158 x 298 | GL45 | 2 |
| | 1519033-FLS | PUREGRIP® Reusable Amber Glass Media Storage Bottles with VersaCap® | 5,000 | 183 x 326 | GL45 | 1 |
| | 1519038-FLS | PUREGRIP® Reusable Amber Glass Media Storage Bottles with VersaCap® | 10,000 | 235 x 408 | GL45 | 1 |
| | 1519040-FLS | PUREGRIP® Reusable Amber Glass Media Storage Bottles with VersaCap® | 20,000 | 296 x 503 | GL45 | 1 |

*EZBio® Cap Assemblies are compatible with the GL45 VersaCap® and can be purchased separately

